


## **Counting chamber, Improved Neubauer**

## **Counting chamber, Improved Neubauer**

With double ruling. Calibratable.

Cell depth 0.1mm. Counting grid with 9 large squares each 1 mm². The 4 large squares in the corners for counting leucocytes are each subdivided into 16 squares with a side length of 0.25mm. The large square in the middle is subdivided into 25 squares each with a side length of 0.2mm. Each of these is further subdivided into 16 smaller squares each with a side length of 0.05 mm and an area of 0.0025mm². The 5 group squares marked with an 'E' are used to count thrombocytes and erythrocytes. All group squares have triple boundary lines on each side.

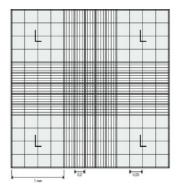



Code	Description	Packaging
LLG07020272	Counting chamber, Improved Neubauer, Type With snap-on clips	1 pz
LLG09161078	Counting chamber, Improved Neubauer, Type Without snap-on clips	1 pz.



CARLO ERBA Reagents S.A.S. Chaussée du Vexin – 27106 Val de Reuil cedex

N° TVA: FR 63391048824 Tél: +33 (0)2 32 09 20 00 www.carloerbareagents.com



CARLO ERBA Reagents operates with a Certified Quality Management System

